



GRUPO DE HABLA ESPAÑOLA Y PORTUGUESA DE LA ISFG

GRUPO DE LÍNGUAS ESPANHOLA E PORTUGUESA DA ISFG

Instituto Nacional de Toxicología
y Ciencias Forenses

MINISTERIO DE LA PRESIDENCIA JUSTICIA Y RELACIONES CON LAS CORTES
SERVICIO DE GARANTÍA DE CALIDAD
DEPARTAMENTO DE MADRID
C/ José Echegaray nº 4 - 28232 Las Rozas de Madrid (Madrid)
Tf.+34 91 768 89 19
e-mail: intcf.eiadm@justicia.es

PROFICIENCY TEST

“ANALYSIS OF DNA POLYMORPHISMS IN BLOODSTAINS AND OTHER BIOLOGICAL SAMPLES”

BASIC LEVEL**EXERCISE EIADN -34 (2026)****DEADLINE: 15/05/2026**

Items sent

2026/Kinship Module

M1 to M3: reference items

2026/Forensic module

M4: forensic unknown item

M5: hair sample

Seal number

Approach:

2026/Kinship Module – Basic level

Practical Kinship study

- **M1, M2, M3:** reference items for genetic profiling.

Theoretical Kinship study

Participants are asked to solve the proposed theoretical study.

2026/Forensic Module – Basic level

Practical Forensic study

- **M4:** forensic item for genetic profiling.
- **M5:** hair for mitochondrial DNA analysis.

◆ Determine the body fluid component or possible components of the item M4.

◆ Could any of the donors from the reference items M1, M2, M3 have contributed to the item M4?

Theoretical Forensic study


Participants are asked to solve the proposed theoretical study.

Methodology to be used

The analyses will be performed by using the markers and methods chosen by the laboratory or those use of routine or that are being under implementation. The items must be processed as real casework and, if possible, as blind samples.

The National Institute of Toxicology and Forensic sciences is the accreditation holder

**Activities marked are not bound to ENAC accreditation*

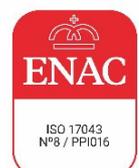

INDEX

	Page
1. Methodology	
1.1. STRs methodology	
1.1.1. Multiplex kits methodology	3
1.1.2. Other methodology for autosomal STR markers and amelogenin	4
1.1.3. Other methodology for Y-STR markers	4
1.1.4. Other methodology for X-STR markers	4
1.2. Mitochondrial DNA methodology	
1.2.1. Amplification parameters	4
1.2.2. Sequencing and editing parameters	5
1.3. Methodology for body fluid identification of item M4	5
1.4. Other considerations regarding methodology different to reported	5
2. Practical studies results	
2.1. STR markers results	
2.1.1. Autosomal STRs and amelogenin	6
2.1.2. Y-STRs	7
2.1.3. X-STRs	8
2.2. Mitochondrial DNA results	8
3. Practical studies conclusions	
3.1. Kinship module	9
3.1.1. Remarks about items M1, M2 and M3	9
3.2. Forensic Module	9
3.2.1. Question 1	9
3.2.2. Question 2	9
3.2.3. Question 3	9
3.2.4. Remarks about items M4 and M5	9
4. Theoretical studies	
4.1. Kinship theoretical study	10
4.2. Forensic theoretical study	12
5. Remarks about this exercise	16
6. Suggestions for subsequent exercises	16
7. Compromises to be met by the participant	16
Date and signature of the person in charge	16
Request for certification	16

1. Methodology *Read carefully the instructions provided before filling in this section***1.1 DNA Extraction, purification/concentration and quantitation**

TABLE 1

Item	Differentia l lysis (Yes/No)	Extraction Purification/ Concentration (Code)	EP00 (Specify)	Quantitation (Code)	C00 (Specify)
M1					
M2					
M3					
M4					
M5					

See Appendix 2026 for codes

1.2 STRs methodology**1.2.1 Multiplex kits methodology**

TABLE 2A (Multiplex kits)

If a kit not included in the table is used, add it in the last rows.

Multiplex	Report 'YES' if used	Detection (Code)	D00 (Specify)
FFFL (Promega)			
PowerPlex 16/16 HS (Promega)			
PowerPlex ESI 16 (Promega)			
PowerPlex ESX 16 (Promega)			
PowerPlex ESI 17 (Promega)			
PowerPlex ESX 17 (Promega)			
PowerPlex 18D (Promega)			
Profiler Plus (AB)			
SGM Plus (AB)			
Identifiler (AB)			
Identifiler Plus (AB)			
Identifiler Direct (AB)			
NGM (AB)			
NGM SElect (AB)			
MiniFiler (AB)			
Investigator ESSplex (Qiagen)			
Investigator ESSplex SE (Qiagen)			
Investigator IDplex (Qiagen)			
YFiler (AB)			
PowerPlex Y (Promega)			
Argus X-8 (Biotype)			
Investigator Argus X-12 (Qiagen)			
XSTR-Decaplex GHEP (Gusmão)			
PowerPlex CS7 (Promega)			
Profiler (AB)			
Investigator Argus Y-12 (Qiagen)			
SEfiler (AB)			
PowerPlex 23Y(Promega)			
PowerPlex Fusion System (Promega)			
Global Filer (AB)			
PowerPlex 21 (Promega)			
Investigador 24plex QS (Qiagen)			

Multiplex	Report 'YES' if used	Detection (Code)	D00 (Specify)
PowerPlex Fusion 6C System (Promega)			
Verifiler (AB)			
Yfiler Plus (AB)			
Investigator ESSplex plus (Qiagen)			
Investigator ESSplex Plus (Qiagen)SE			
Investigator IDplex Plus (Qiagen)			
Investigator HDplex (Qiagen)			
Investigator Argus X-12 QS (Qiagen)			

See Appendix 2026 for codes

1.2.2 Other methodology for autosomal STR markers and amelogenin

TABLE 2B

Fill in ONLY in case kits multiplex are not used or additional STR markers are utilized. Indicate the number of markers, the primers and the methodology used.

Number of markers	Primer/Ladder (Code)	PL00 (Specify)	Detection (Code)	D00 (Specify)

See Appendix 2026 for codes

1.2.3 Other methodology for Y-STR markers

TABLE 2C

Fill in ONLY in case kits multiplex are not used or additional Y-STR markers are utilized. Indicate the number of markers, the primers and the methodology used.

Number of markers	Primer/Ladder (Code)	PL00 (Specify)	Detection (Code)	D00 (Specify)

See Appendix 2026 for codes

1.2.4 Other methodology for X-STR markers

TABLE 2D

Fill in ONLY in case kits multiplex are not used or additional X-STR markers are utilized. Indicate the number of markers, the primers and the methodology used.

Number of markers	Primer/Ladder (Code)	PL00 (Specify)	Detection (Code)	D00 (Specify)

See Appendix 2026 for codes

1.3 Mitochondrial DNA methodology

1.3.1 Amplification parameters

TABLE 3

Report each primer set in one single box and name them according to the strand (L or H) and 3' nucleotide position (Ex. L15997/H00619)

Primers sets for amplification					Nr of cycles
Item	Forward/reverse	Forward/reverse	Forward/reverse	Forward/reverse	
M1-M3					
M4					
M5					

2. Practical studies results:

Read carefully the instructions provided in order to fill in the results tables and the rules of participation in order to know the establishment of assigned values and the evaluation of results
<https://ghep-isfg.org/en/proficiency/participation/>

2.1 STRs Results

ALL PARTICIPANTS OF THE FORENSIC MODULE, MUST COMPLETE COMPULSORY THE COLUMN OF TOTAL ALLELES DETECTED REGARDLESS THE EXTRACTION SYSTEM USED. The 1st and 2nd fraction columns are additional and optional, in case the laboratory have performed differential lysis and want to reflect its result.

2.1.1 Autosomal STRs and amelogenin

TABLE 6A

KINSHIP MODULE				FORENSIC MODULE		
MARKER	M1	M2	M3	M4		
				Total of alleles detected Ex:9-11-15-17	1 ST fraction Ex:9-17	2 nd fraction Ex: 11-15
AMEL						
D8S1179						
D21S11						
D7S820						
CSF1PO						
D3S1358						
TH01						
D13S317						
D16S539						
D2S1338						
D19S433						
vWA						
TPOX						
D18S51						
D5S818						
FGA						
Penta D						
Penta E						
D10S1248						
D22S1045						
D2S441						
D1S1656						
D12S391						
SE33						
FES/FPS						
F13A01						
F13B						
LPL						
Penta C						
D6S1043						

2.1.2 Y-STRs

TABLE 6B

KINSHIP MODULE				FORENSIC MODULE		
MARKER	M1	M2	M3	M4		
				Total of alleles detected Ex: 13-15	1 ST fraction Ex: 15	2 nd fraction Ex: 13
DYS456						
DYS389 I						
DYS390						
DYS389 II						
DYS458						
DYS19						
DYS385						
DYS393						
DYS391						
DYS439 (GATA A4)						
DYS635 (GATA C4)						
DYS392						
GATAH4						
DYS437						
DYS438						
DYS448						
DYS460 (GATA A7.1)						
DYS461 (GATA A7.2)						
GATAA10						
DYS388						
DYS576						
DYS481						
DYS549						
DYS533						
DYS570						
DYS643						
DYS627						
DYS518						
DYS449						
DYF387S1						

2.1.3 X-STRs

TABLE 6C

KINSHIP MODULE				FORENSIC MODULE		
MARKER	M1	M2	M3	M4		
				Total of alleles detected Ex: 12-15-17-20	1 ST fraction Ex: 12-15	2 nd fraction Ex: 17-20
HPRTB						
DXS8378						
DXS9898						
DXS7133						
GATA31E08						
GATA172D05						
DXS7423						
DXS6809						
DXS7132						
DXS9902						
DXS6789						
DXS10103						
DXS10134						
DXS10074						
DXS10101						
DXS10135						
DXS10146						
DXS10079						
DXS10148						

2.2 Mitochondrial DNA results

In Table 7A, report the initial and final positions of the edited regions and in Table 7B report the haplotypes in the order requested in the instructions

TABLE 7A

ITEMS		EDITED REGIONS
KINSHIP MODULE		
M1		
M2		
M3		
FORENSIC MODULE		
M4	1 st fraction	
	2 nd fraction	
Hair M5		

TABLE 7B

ITEMS		HAPLOTYPE
KINSHIP MODULE		
M1		
M2		
M3		
FORENSIC MODULE		
M4	1 st fraction	
	2 nd fraction	
Hair M5		

3. Practical Studies Conclusions

3.1 Kinship Module

3.1.1 *Remarks about items M1, M2 and M3

Indicate any comments or remarks, you consider, about the analyzed items. Please, remember that only the genetic profiling of the reference items M1 to M3 is required; it is not necessary to investigate a genetic relationship among them.

3.2 Forensic Module

3.2.1 Determine the body fluid component or possible components of the item M4.

Components (mark with an X the component/s detected)

	Blood	Semen	Saliva
M4	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

3.2.2 Indicate the minimum number of contributors detected in the item M4.

	1	2	3
M4	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

3.2.3 Could any of the donors from the reference items M1, M2, M3 have contributed to the item M4?

	M1	M2	M3
M4	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

3.2.4 *Remarks about items M4 and M5.

4. Theoretical studies

Read carefully the instructions provided in order to fill in the results tables and the rules of participation in order to know the establishment of assigned values and the evaluation of results
<https://ghep-isfg.org/en/proficiency/participation/>

In order to solve the theoretical studies (kinship and forensic) it is assumed that:

- the population is in Hardy-Weinberg equilibrium and that no correction is made due to population substructure ($\theta=0$).
- silent alleles rate and mutation rate are 0.
- drop in, drop out correction=0

Calculations have to be made by using the “2026 Alleles Frequencies” table provided.

4.1 Kinship theoretical study

4.1.1 Approach

Albert Jr. claims to be the son of a well-known businessman, Mr. Albert, who has passed away. His mother, Ms. Pilar, worked for many years as the Company's secretary, and she was dismissed when she became pregnant by the businessman. Nevertheless, the businessman continued to pay child support until Alberto Jr. reached the legal age.

The judge finds enough evidences to accept the paternity suit and orders a DNA test. The alleged father's body was cremated, but a Formaline-Fixed Paraffin-Embedded (FFPE) tissue of a coronary artery biopsy, performed at the hospital before his death, is available. A genetic profile is obtained from the FFPE tissue.

- The genetic profile of Ms Pilar (Albert Jr's mother) s is available. There is no doubt that Ms Pilar is Albert Jr's biological mother.
- The genetic profile of Albert Jr is available
- A partial genetic profile obtained from the FFPE tissue is available (Mr Albert, alleged father)

Markers	Ms. Pilar	Albert Jr.	Alleged father (FFPE)
D3S1358	15	14-15	14-16
D1S1656	12-15	12	12-18.3
D2S441	11-13	11-14	10-14
D10S1248	15	15	
D13S317	8-14	8-9	
D16S539	11-12	11-12	9-11
D18S51	13-17	14-17	12-14
D2S1338	19-23	17-19	
Penta D	10-15	10-15	12-15
TH01	7-9	7-8	8-9.3
VWA	15-18	18-19	17-19
D21S11	31.2-32	29-32	28-29
D8S1179	13-14	12-13	12-14
D19S433	13-16	13-15.2	14-15.2
D22S1045	14-15	14-16	15-16
FGA	19-26	21-26	
AMELOGENINE	X	X-Y	X-Y

4.1.2. Paternity index

Calculate the paternity index taking into account the following hypotheses:

H0	Mr. Alberto is the father of Alberto Jr., considering that Mrs. Pilar is the biological mother.
H1	Another man randomly selected from the population and not genetically related to Mr. Alberto is the biological father of Alberto Jr., considering that Mrs. Pilar is his biological mother.

Report the partial paternity indexes and the total PI in **Table 8**.

Use scientific notation (Excel format) and rounding off to 4 decimals places Ex. 1,2346E-01

TABLE 8

Markers	PI
D3S1358	
D1S1656	
D2S441	
D10S1248	
D13S317	
D16S539	
D18S51	
D2S1338	
Penta D	
TH01	
VWA	
D21S11	
D8S1179	
D19S433	
D22S1045	
FGA	
TOTAL PI	

4.1.3 Software/s used to carry out the statistical calculations.

Program	version	Remarks (other software, comments, etc)
Familias		
DNA view		
PatPCR		
BDGen		
PatCan		
Genética Forense Final		
Home-made Software		
Others ¹		

¹If your software is not displayed in the table, chose "others" and specify it in the cell Remarks

4.1.4 Hand-made calculations. Formulas used

In the case, your laboratory performs all calculations by hand, then report the used formulas in **Table 9**

TABLE 9

Markers	PI
D3S1358	
D1S1656	
D2S441	
D10S1248	
D13S317	
D16S539	
D18S51	
D2S1338	
Penta D	
TH01	
VWA	
D21S11	
D8S1179	
D19S433	
D22S1045	
FGA	
TOTAL PI	

4.1.5 *Conclusions and remarks about the kinship theoretical study.

4.2 Forensic theoretical study

4.2.1 Approach

A young woman with disabilities reported a sexual assault perpetrated by several individuals. he states that her boyfriend (with whom she had never had sexual relations), his brother and a cousin of both of them participated in the assault.

Following an examination at the hospital, a vaginal wash is performed, after analyzing a mixed DNA profile is obtained.

The boyfriend (N), the boyfriend's brother (H), and their cousin (P) are arrested. Reference DNA samples are taken. **The reference profile of the victim is available.** The corresponding genetic profiles are obtained from the reference samples.

Determine if any of the three males are compatible with the mixture and, based on that, calculate the LR.

Marcadores	Vaginal wash Mixed genetic profile	Victim (V)	Boyfriend (N)	Boyfriend's brother (H)	Their cousin (P)
D3S1358	15-16	15-16	15-16	15-16	15-16
D1S1656	11-15-16-18.3	15-16	12-18.3	11-18.3	16-17.3
D2S441	11-14-15	14-15	10-11	11	10-11
D10S1248	13-15	13	13-15	15	13-15
D13S317	9-12-13	12-13	8	9-13	11-13
Penta E	5-12-13-14	5-14	12-13	12-13	13-17
D16S539	9-11-13	9-11	9-12	9-13	9-13
D18S51	14-15-16	14	15	15-16	14
D2S1338	17-19-20-24	19-24	17-20	17-20	20-24
CSF1PO	10-11-12	11	9-10	10-12	11
Penta D	9-13-15	9-13	9-10	9-15	9-16
TH01	6-7-9.3	6-9.3	6-7	6-7	6-9.3
vWA	17-18-19	18-19	18	17-19	15-17
D21S11	29-30.2-32.2-33.2	29-32.2	30-33.2	30.2-33.2	28-31.2
D7S820	8-10-11-12	8-10	11-12	11-12	8-10
D5S818	9-10-11	10	9-12	9-11	7-14
TPOX	8-10-12	10-12	8-11	8-12	8-9
D8S1179	12-13-15-16	12-13	12-14	15-16	14-15
D12S391	15-17.3-18	15-18	18-20	17.3-18	18
D19S433	13-14-14.2-16.2	14.2-16.2	14	13-14	13-14
SE33	14-15-16-19	15-16	11-17	14-19	15.2-20
D22S1045	14-15-16	14-15	15-16	15-16	16
FGA	19-21-22-23	19-22	21-23	21-23	21
Amelogenine	X-Y	X	X-Y	X-Y	X-Y

4.2.2. LR value

Report the partial Likelihood Ratio (LR) values, as well as the total LR in **Table 10**, according to the following hypothesis:

H0	The victim and the person involved (N, H or P) have contributed to the mixed genetic profile obtained in the vaginal wash.
H1	The victim and an unknown person taken at random from the population and not genetically related to the previous ones have contributed to the mixed genetic profile.

Indicate if one of these males is present in the mixed genetic profile*

Boyfriend (N)		Boyfriend's brother (H)		Their cousin (P)	
---------------	--	-------------------------	--	------------------	--

TABLE 10

Use scientific notation (Excel format) and rounding off to 4 decimals places . Ex. 1,2346E-01

Marcadores	LR
D3S1358	
D1S1656	
D2S441	
D10S1248	
D13S317	
Penta E	
D16S539	
D18S51	
D2S1338	
CSF1PO	
Penta D	
TH01	
vWA	
D21S11	
D7S820	
D5S818	
TPOX	
D8S1179	
D12S391	
D19S433	
SE33	
D22S1045	
FGA	
Total LR	

4.2.3 Software/s used to carry out the statistical calculations.

Program	version	Remarks (other software, comments, etc)
LRmix Studio		
LR mezcla v inteligente		
EuroForMix		
DNAMix		
Genética Forense Final		
Home-made Software		
DNA View		
Others ²		

²If your software is not displayed in the table, choose "others" and specify it in the cell "Remarks".

4.2.4 Handmade calculations. Formulas used.

In case of only handmade calculations, then report the used formulas in **Table 11**.

TABLE 11

Marcadores	LR
D3S1358	
D1S1656	
D2S441	
D10S1248	
D13S317	
Penta E	
D16S539	
D18S51	
D2S1338	
CSF1PO	
Penta D	
TH01	
vWA	
D21S11	
D7S820	
D5S818	
TPOX	
D8S1179	
D12S391	
D19S433	
SE33	
D22S1045	
FGA	
Total LR	

4.2.5 *Conclusions

Issue a conclusion regarding the results obtained.

--

5. Remarks about this exercise

6. Suggestions for subsequent exercises

7. Compromises to be met by the participant

The analyses, both, the generated results and their statistical evaluation have been performed in the facilities of the participating laboratory and by its own staff, following working protocols used in routine casework together with safety precautions. **In accordance with the donors' consent, these items will be processed anonymously for the Intercomparison Exercise INTCFM/GHEP-ISFG. Additionally they could be used as a reference material and/or quality control for the laboratory either using the techniques required in the Exercise or other forensic techniques but always for the purpose of human identification, analyzing non coding regions or regions that would not provided sensitive information about the donor: illnesses, pathologies or other genetic information which could infringe his/her privacy.**

Name of the person in charge

Date and signature

WOULD YOU LIKE TO RECEIVE A CERTIFICATE OF PARTICIPATION?

Kinship Module (Basic level) (Yes/No)	Forensic Module (Basic level) (Yes/No)

CHOOSE THE LANGUAGE OF THE CERTIFICATE

SPANISH
 ENGLISH
 BOTH

Note.- In order to receive the certificate of participation you must return this form duly signed.